



# The New Algorithm for Knee Cartilage Treatment



Biological Grafting

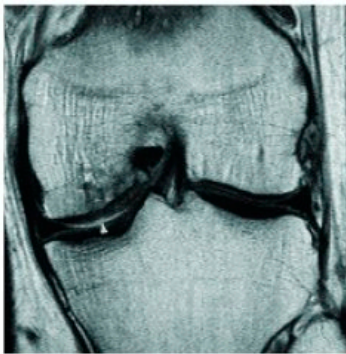

AKR Knee Resurfacing

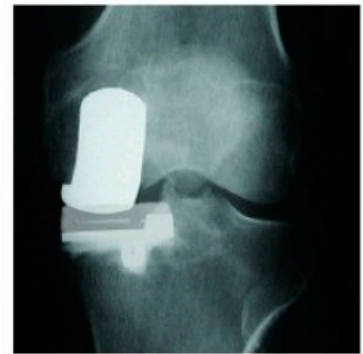

Unicompartmental Knee Replacement



### Femoral Component

- Cobalt chrome articulating surface
- Titanium plasma spray undercoating
- Bead-blasted titanium taper post

### Tibial Component

- Inlay cemented UHMWPE

The AKR surgical technique utilizes a cannulated system similar to those used in knee ligament repair.

### Tibial



Tibial Guide



Retrograde Socket Preparation



Sizing Trial




Tibial Implant


### Femoral



Femoral Guide

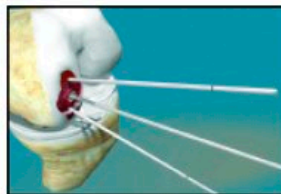

Set Guide Wires

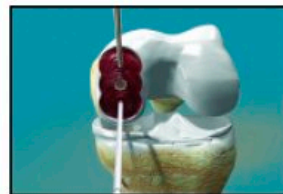

Triple Ream Femoral Sockets



Femoral Implant

PN 2001-4002 REV C

28 forge parkway, franklin, ma 02038 tel +1 508 520 3003 fax +1 508 528 4604 web arthrosurface.com  
 This product is covered by one or more of US Patent Nos. 6,520,964, 6,610,067, 6,679,917 and other patents pending.  
 UniCAP™ is a trademark of Arthrosurface®, Inc.